


# Open Rocker Switch

## QUICK GUIDE

### Introduction

The Open Rocker Switch is an accessible switch with two buttons. The switch is activated by pressing down on one side or the other like a see-saw. This device is suitable for anyone and can be connected to any device with a 3.5 mm interface.

### Features

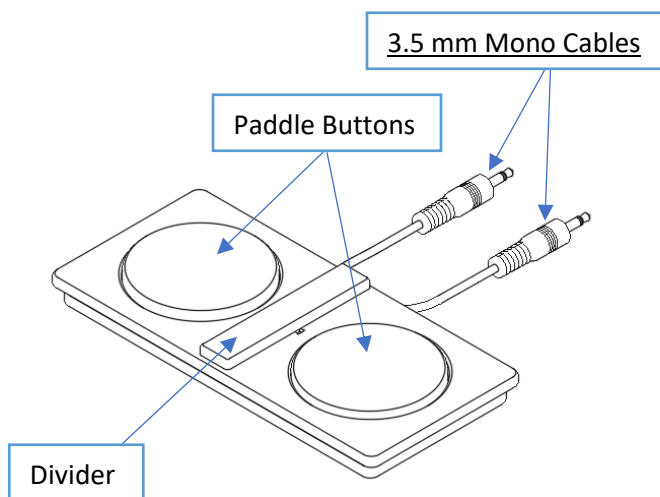

### Usage

Connect the cables to the assistive device.

Action	Action Time (seconds)	Function
Press Button Left	< 1	Depends on interfaced device
Press Button Right	< 1	Depends on interfaced device

### Compatibility

Compatible with any device that accepts a 3.5 mm mono cable. The device can be connected to two separate devices, one device with two mono inputs, or one device with a stereo input with an appropriate mono to stereo adapter.

### Specifications

<b>Activation Force</b>	1.27 ± 0.49 N 130 ± 50 grams-force
<b>Activation Distance</b>	~ 0.3 mm
<b>Dimensions (Length x Width x Height)</b>	170 x 80 x 18.35 mm
<b>Weight</b>	~100 grams (Depends on print material)
<b>Switch Electrical Rating</b>	50 mA at 24VDC

### Cleaning

- Disassemble all pieces
- Wash paddle buttons, divider, and top piece in any cleaning solution
- Wipe clean bottom piece with a towel or paper towel and leave to dry completely to ensure that the buttons will not short circuit